

WANROOEMED

Your Ideal Operating Table Supplier

Well Sold For 40 Years Worldwide
Meet The Basic Demands Of General Surgery

AOT700 Electro-hydraulic operating table

With reference to the advantages of congener products from home and aboard,can meet the demands of different position on operation.
24V battery that meets the operating demands for one week,better ensure the table work under safe voltage.
Electromotion facilitates the various adjustment of position,intellectualized reposition enables the table a quick return to horizontal position.Automatically lockable remote control in case of wrong operation.
Large-scale elevation makes it easier for C-arm fluoroscope examination.
With built-in kidney bridge,suits for kidney operation.
Optional traction for orthopedics,head rest for Aural surgery ,extend functions of the table .The cover all adopted premium 304 stainless steel,the main body used cast aluminum,anti-rust and corrosion the table top adopted phenol aldehyde photograph board which is suitable for X-ray.
Key parts such as electrical pump and electromagnetic valve,adopted imported products from Taiwan.



Your Ideal Operating Table Supplier

Well Sold For 40 Years Worldwide
Meet The Basic Demands Of General Surgery



江苏如翼医疗科技有限公司
JIANGSU ROOE MEDICAL TECHNOLOGY CO.,LTD.
Add: Room1311,1313, Guotai New World Plaza, No19,
East Renmin Road,Zhangjiagang, Jiangsu, China.
Tel:+86 512 58810663
Whatsapp/Mobile:+ 86-18962287398;15962399164
Email:info@wanrooemed.com
Website:www.rooemed.com www.wanrooemed.com

WANROOEMED

Your Ideal Operating Table Supplier

Well Sold For 40 Years Worldwide
Meet The Basic Demands Of General Surgery

AOT700

Electro-hydraulic
operating table



JIANGSU ROOE MEDICAL TECHNOLOGY CO.,LTD.

WANROEMED

Your Ideal Operating Table Supplier

Well Sold For 40 Years Worldwide
Meet The Basic Demands Of General Surgery



PARAMETER A0T700 Electro-hydraulic operating table

Movements		Model: A0T700
Tabletop Reverse Trendelenburg		≥28°
Tabletop Trendelenburg		≥28°
Tabletop Lateral Tilt (Left And Right)		≥18°
Head Plate	Up	≥60°
	Down	≥85°
Back Plate	Up	≥85°
	Down	≥40°
Leg Plate	Up	≥15°
	Down	≥90°
	Outward	≥90°
Horizontal Sliding		≥340mm
Kidney Bridge Up		100mm

WANROEMED


Your Ideal Operating Table Supplier

Well Sold For 40 Years Worldwide
Meet The Basic Demands Of General Surgery

Model	Site	Length	Width	Tabletop Lowest Height	Elevation range
A0T700		≥1960mm	500mm	≤690mm	400mm

STANDARD ACCESSORIES A0T700 Electro-hydraulic operating table

Anesthetic Screen	Shoulder Support	Arm Support	Body Support
1piece	1pair	1pair	1pair
Leg Support	Fixing Clamp	Foot Plate	Remote Control
1pair	1pair	1pair	1piece
Mattress	Kidney Bridge Handle	Power Line	Hydraulic Oil
1 set	1piece	1piece	1oilcan

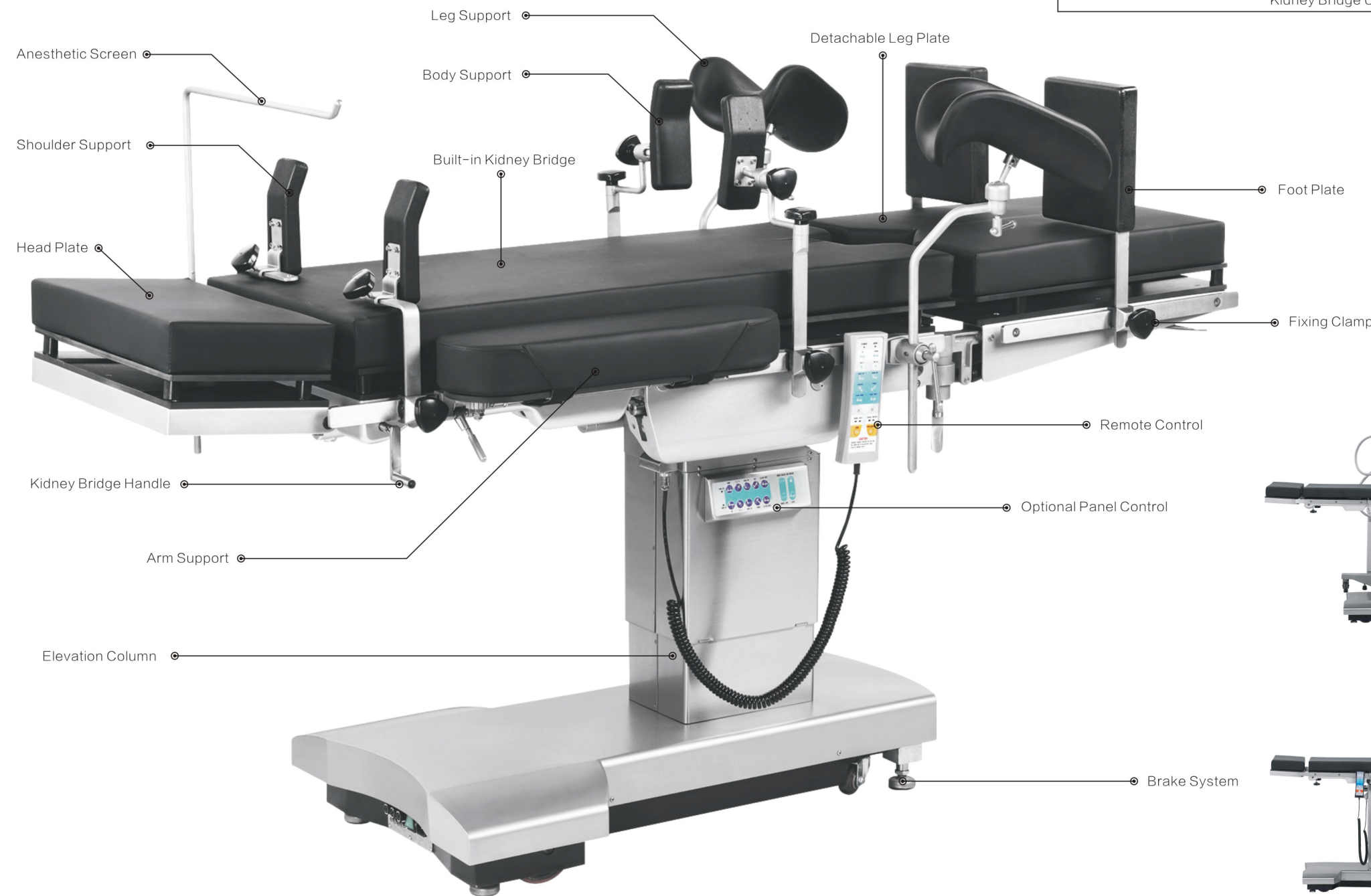

Your Ideal Operating Table Supplier

Well Sold For 40 Years Worldwide
Meet The Basic Demands Of General Surgery

