



#### More classic practical function

- Scanning mode: B、2B、4B、B/M、M
- Full field image with 8 levels depth adjustment
- Dynamic range: 27-90dB, visual adjustment
- Local zooming function
- Image storage: 100 frames permanent storage, without loss of the image when powered off
- Cine loop: 256 frames
- High-speed USB-A port: Supporting Large-capacity USB flash disk and mouse
- High-speed USB-B port: for image transfer to PC (Optional Case Report Processing Software)


#### Professional measuring software for veterinary

- Common measurement: distance, circumference, area, volume, etc.
- OB measurement: Bovine, Equine, Ovine, Canine, Feline, Goat, Llama, gestational weeks and expected date.
- Cardio measurement: depth, slope, heart rate, cycle

#### Configuration

- Standard configuration: (1) Main unit : 1pc (2) 6.5MHz rectal vet probe (multi frequency, without handle) : 1pc (3) Internal lithium battery: 2 pcs (4) Charger: 1pc (5) Belt: 1pc (6) Bracket: 1pc (7) Leather case : 1pc (8) Handle freeze switch
- Optional Configuration: 4.0MHz multi-frequency rectal convex Probe、 6.5MHz multi-frequency linear rectal probe (handle)、 3.5MHz multi-frequency convex probe、 7.5MHz high multi-frequency liner probe、 5.0MHz multi-frequency micro convex probe、 car charger、 goggle, etc

#### Image quality

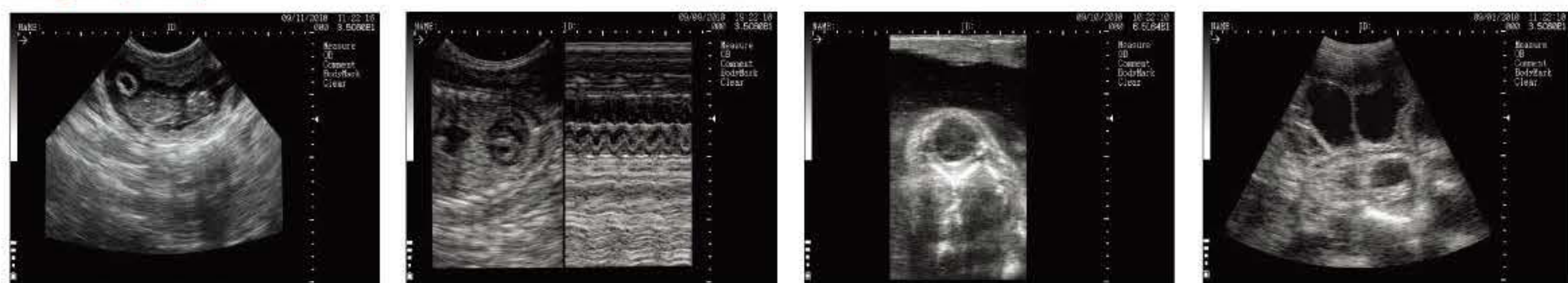

Gravidity of dog which is 30 days

Heart of dog which is 50 days

Eyeball of cow

Gravidity of pig which is 25 days

#### Optional probe



4.0MHz multi-frequency rectal convex Probe



6.5MHz multi-frequency linear rectal probe (handle)



3.5MHz multi-frequency convex probe



7.5MHz high multi-frequency liner probe



5.0MHz multi-frequency micro convex probe



万容机械有限公司

WANROOE MACHINERY CO., LTD.

江苏省张家港市人民东路19号国泰新天地广场1311,1313室

Room1311,1313, Guotai New World Plaza, No19,

East Renmin Road,Zhangjiagang, Jiangsu, China.

Tel:+86 512 58810663

Fax:+ 86 512 58810663

Email:doris@wanrooemed.com

DISTRIBUTOR

THE PRODUCT TAKES REAL OBJECT AS QUASI,IF ANY FLUCTUATION HAPPENS,FORGIVING NO ANOTHER NOTICE IN THE LINE.

# USKX5200V

Full-digital B mode Veterinary Ultrasound Scanner



## Features



Two large-capacity lithium batteries, about 360 minutes work.



8 kinds of internal Pseudo color display



5.7 Inch high resolution colour TFT LCD monitor forming a clear image under the sunshine for special treatment.



1.1Kg

Weight of main unit is 1.1Kg, small and light



Backpack integration design, convenient for veterinary and farm use



The main unit with the splash-proof, dust-proof function, easy to clean.



#### Main unit

Smart design and internal battery of USKX5200 make the diagnosis more convenient and bring your new mobile using experience.

Create the health with science and technology