


Advantages:

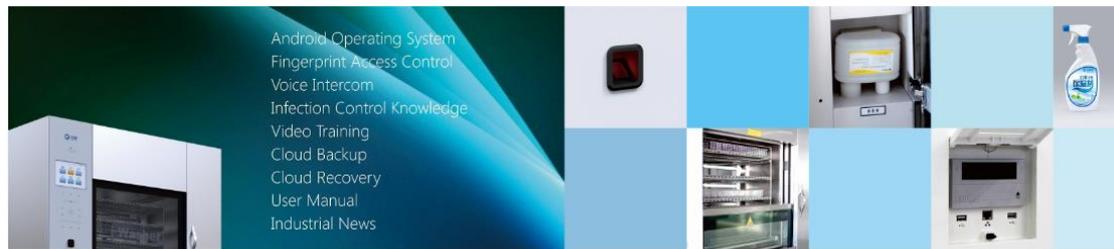

- The first manufacturer in the world to implement Android operating system to washer machine, which indicates that infection control equipment break into informational and intelligent age.
- The first manufacturer in the world to apply fingerprint access control to washer machine, which is easier and faster for operators to enter the interface, and enhances the traceability.
- The first washer machine in the world to be equipped with cloud service enabling real-time monitoring anytime and anywhere.
- The first washer machine in the world to have the function of text message notification (SMS), informing the operator of running status, the cycle end, equipment faults and maintenance, etc.
- The first and only one in the world to play training video, online update of industrial news, infection control knowledge, and user manual.
- The first washer machine in the world with the function of voice intercom enabling real-time communication between decontamination area and clean area, which provides an easy way to communicate for CSSD staff.

Quality:

- CE
- FDA
- FSC
- ISO9001
- ISO13485

Patents:

- Patent of high-speed automatic sliding door
- Patent of high-speed elastic connector
- Patent of water cooling and drying device for washer disinfector
- Patent of washer disinfector with traceability
- Patent of washer disinfector with fingerprint access control
- Patent of washer disinfector with SMS
- Patent of washer disinfector with voice intercom system

- **Smart operating system:**

to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663



WANROOEMED

The equipment implements Android operating system, enabling the real-time monitoring of machine. The running status and record can be delivered to traceability system by wire or wireless connector. In this way, the washer disinfector is monitored real time.

- **Finger print access control:**

Users can enter operating interface by finger print ID or password. By using finger print ID technology, the equipment gets physiological feature record of operators, so it can trace the responsibility of operators.

- **Video training:**

User can receive operation training and infection control knowledge, etc from the machine.

- **Text message notification(SMS):**

The equipment operating staff can get the information about running status of machine by text message. (including running information and fault information)

- **Various connectivity to network:**

Provide RJ 45 wire connector, WIFI module, 3G network module and USD data interface.

- **Cloud service:**

After connection with network, the equipment can backup the data into cloud on Internet. Once there is any fault, the machine, after repaired, can synchronize data from cloud.

Any operators who are authorized can monitor the running of machine anytime and anywhere.

Main functions of cloud: cloud back up, cloud recovery, cloud curve.

- **Online service:**

Provide the latest news in medical field;

Launch the latest equipment information;

Update the newest industry knowledge;

Provide operating manual and maintenance knowledge;

Online service is available anytime.

- **Voice intercom system:**

BAIX voice intercom system is recognized by China patent department, which enables the real-time communication between decontamination area and clean area.

The function of voice filter can identify voice between staff and environment, so that the high quality voice is guaranteed.

- **Double synchronized high definition LCD:**

10.4 inch LCD with image resolution 800*600; user friendly interface.


- **Standard spare parts:**

All major parts sourced from professional manufacturers, parts life is extended due to short cycle time; high efficiency and longer using life.

to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663

Full automatic washer disinfector QPQ-500/500

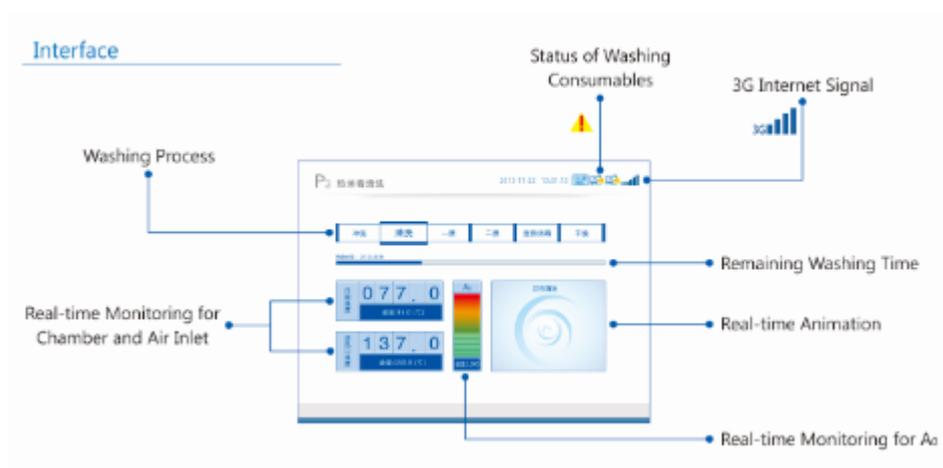

China CDC(Centre for disease control) Test

Washing effect:

China CDC test QPQ 500 washer and disinfector by washing common surgical instruments, dental drill, anesthesia and breathing tubes, humidification bottles, and baby

to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663

feeding bottle. The result of absorbance values of residual protein pollutant is less than 0.020, and approval rate is 100%.

Disinfection programme:

China CDC test QPQ 500 washer and disinfector according to sterilization technical standards and programmes, including low temperature disinfection, high temperature disinfection and chemical disinfection, are approval.

Technical specification:

Model	Power supply	Temperature of heat sterilization	Capacity	Water consumption/cycle	Preset programmes	Overall dimensions	Power
QPQ500	AC380V 50Hz	0-98℃	500L	30L	6	1870*1340*1026(mm)	35KW
QPQ550	AC380V 50Hz	0-98℃	550L	30L	6	1870*1340*1026(mm)	35KW

Racks for Full Automatic Washer Machine



to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663

Full automatic washer machine QPQ-320/360



China CDC(Centre for disease control) Test

Washing effect

China CDC test QPQ 320/360 washer and disinfecter by washing common operation instruments,dental drill,anesthesia and breathing tubes,humidification bottles,and baby feeding bottle.The result of absorbance values of residual protein pollutant is less than 0.020,and approval rate is 100%.

Disinfection programme

China CDC test QPQ 320/360 washer and disinfecter according to sterilization technical standards and programmes,including low temperature disinfection,high temperature disinfection and chemical disinfection,are approval.

Certificates and patents:

Patent of high-speed elastic connector

Patent of water cooling and drying device for washer disinfecter

Patent of washer disinfecter with traceability

Certificates of ISO13485,ISO9001,and CE,etc.


to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663

Technical specification:

Model	Power supply	Temperature of heat sterilization	Capacity	Water consumption/cycle	Preset programmes	Overall dimensions	Power	Basket
QPQ 360	AC380V 50Hz	20-98°C	360L	15L	6	1865*700*700(mm)	17Kw	10
QPQ 320	AC380V 50Hz	20-98°C	320L	15L	6	1860*700*700(mm)	17Kw	8

Operation Instructions

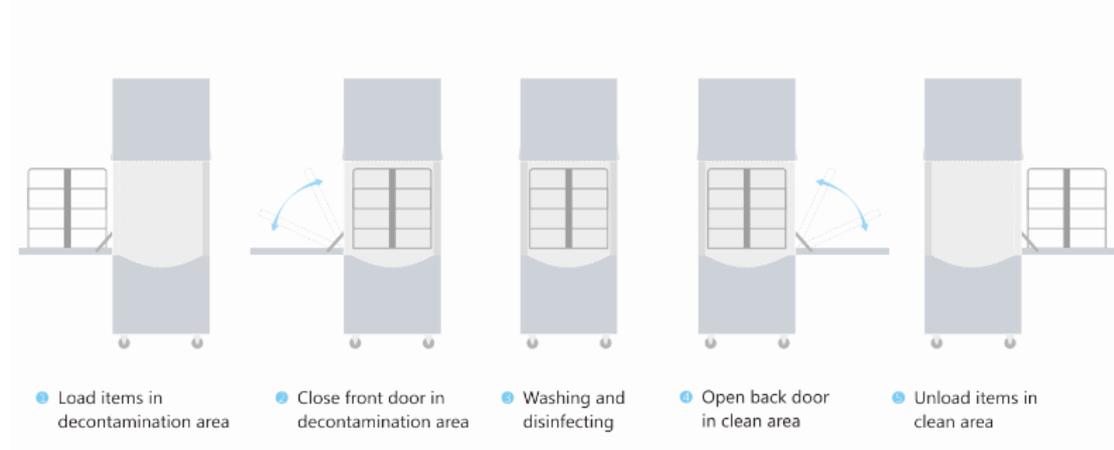

to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663

Full automatic washer disinfector QPQ-200



Product feature:

Humanized control system and easy operation;

Single door,small footprint;

Shorter cycle time and higher efficiency.

Range of application:

Surgical instruments,dental equipment,and endoscooper,etc.

Larger capacity:

Instruments of 2 operation and 24 dental drills.

to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663

Full automatic washer disinfecter

QPQ-60 (Bench top)



Product feather:

Bench top,space-saving

Small size and it can be easily moved

Shorter cycle time and higher efficiency

Range of application:

Surgical instruments,dental equipment,and endoscooper.etc.

Capacity:

Instruments from 1 operation and 12 dental drills

Technical specification:

Model	Power supply	Capacity	Preset programmes	Overall dimension
QPQ-60	AC220V 50Hz	60L	P1:standard washing(28min).P2:quick washing(11min)	570*470*500(mm)

to order

All specifications, pictures are subject to change without earlier notification.
www.wanrooemed.com TEL: 0086-512-58810860 FAX: 086-512-58810663